

Charging terminal

15

case of missing accessory sets, please contact customer service

Care for the top and back of the head.

The correct way to use the brush

The entire brush touches the head

As shown in the figure, move from the top of the head to the hairline in a sequence of 5 seconds each time.

* It is difficult for the scalp to feel the stimulation of EMS, but it is easier to feel

Methods of use

Scalp and Body

Attachment

Needle-shaped

(for scalp use)Attachment%

*Please do not use the needle-shaped attachment on your face or body. It has a strong grip and may pinch.

Face Attachment

Part of the brush is lifted from the scalp

The main machine has waterproof function. It can be used in the bath, but please do not soak it in the bathtub. Please put it into the water carefully. Otherwise, the machine will stop operating and cannot be used.

(right side part)

Attachments


* If using on the face, please be sure to remove makeup.


* If the skin is too dry it will make it difficult to feel the EMS stimulation To achieve the actual physical stimulation of


* The charger, USB cable and AC adapter do not have waterproof function. Please keep it in a dry and waterless place.

*Do not use it in a damp status or charge it in the bathroom.

If using on the face, please be sure to use low-speed model and do not press hard.
Please notice that when using on the head and neck, the stimulation of EMS will be felt stronger than usual

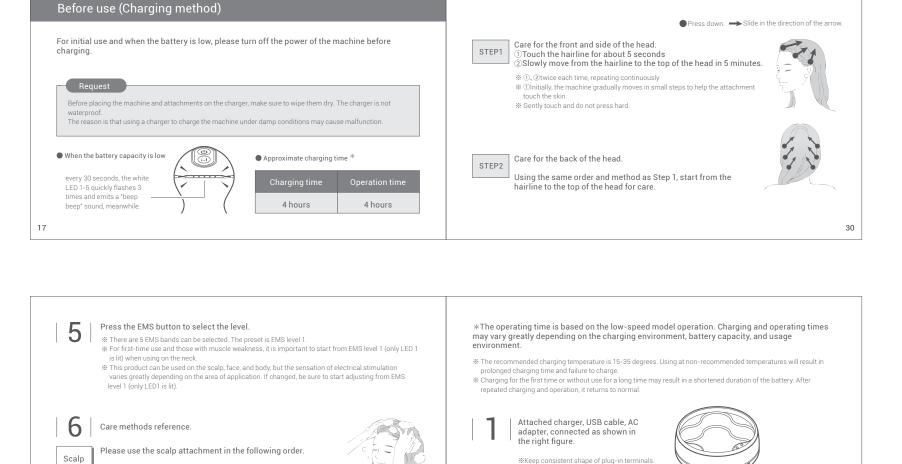

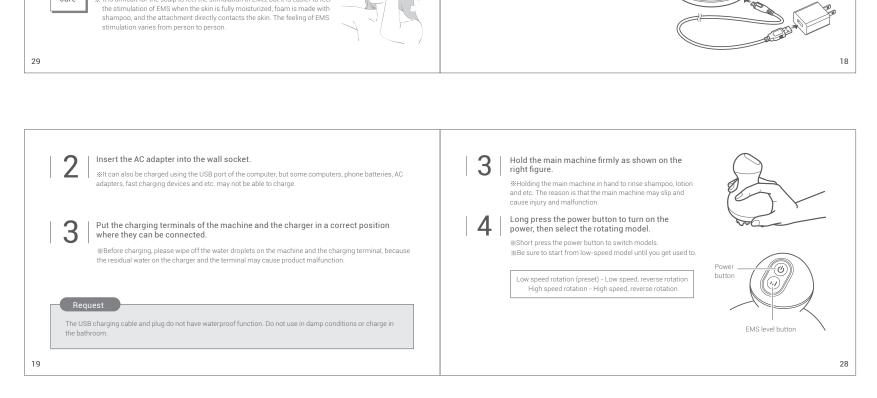

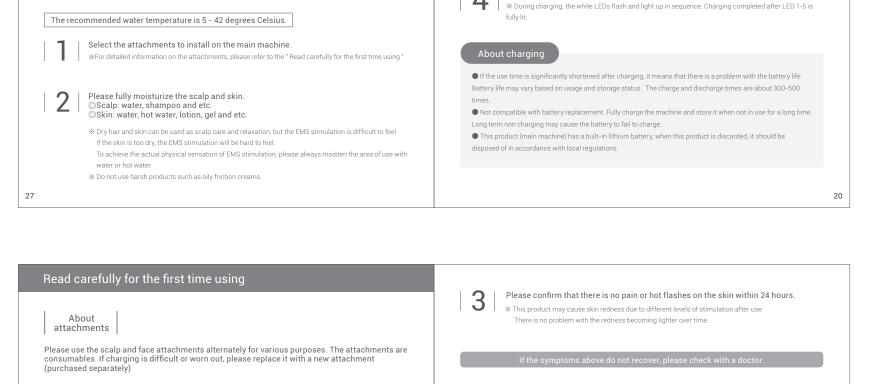

After charging, unplug the AC adapter from the socket.

