

Instruction Manual

Retailer: Sotaro Medical Co., Ltd.
 3-1-10-5F Hakatanishinmachi, Hakata-ku, Fukuoka-shi, Fukuoka-ken
<https://mytrex.jp>


Designing: Japan, Made in China
 ©2022 MYTREX MEDICAL CO., LTD.

MYTREX PROVE
 Full Body Lifting All-in-One Beauty Device


1. Warranty Certificate Attached
 Thank you for purchasing our products. Please read this instruction manual carefully before using to ensure safe and correct usage.
 This user manual comes with a warranty certificate. Please keep it properly after reading.

For users, we recommend watching the instructional video to confirm proper usage of the device.
 Watching the video while using the device can help you achieve more effective results.
 To ensure the video, simply plug the USB cable into...

If you have any questions about the product, please contact customer service.
MYTREX Customer Service
 8:00am to 5:00pm JST, 7/24




Use the Sleep Study device per the targeted stimulation of the eye and mouth muscle.



Use in the order of facial muscles (see page 137-38) especially for targeted care of the eye and mouth area.

Product name
MFTX2 exclusive eye accessory

https://myfreespace.jp/en/
For more information, please visit our official website.



Safety Precautions

Please use the device according to the established procedures and do not use it in an improper way. If any abnormal external event occurs, please stop using the device immediately, and contact the store where you purchased it or our customer service for assistance. Please do not continue to use the device if it malfunctions, as this could have adverse effects on your health and be dangerous.

WARNING

Prohibited

- The following medical electronic equipment should not use the device (inserted part) of the medical electronic device will cause a serious impact on the body.
- Completed medical electronic devices such as pacemakers for the heart.
- Lead-assisted medical electronic equipment such as artificial heart muscle.
- Otherable medical electronic devices such as electrocardiographs.

[illegible]

- The device should not be used in areas other than the designated areas and should not be used in the following areas:
 - Openly flammable or explosive atmospheres
 - Areas where static electricity has been performed
 - Areas with high humidity and high humidity, excessive air pressure, noise, vibration, X-rays, magnetic fields, and so on.
 - Areas with flammable and toxic gases
 - Oil-contaminated tools and work areas
 - Other heat sources (e.g., other electrical tools, etc.) and areas around the wires
 - Other areas around the heat
 - Other obstacles and areas requiring installation
- Devices that have been installed in the same way
 - Devices with integrated heat and oil-proof

WARNING

Do not use the device and wear accessories in an unstable situation, and prevent them from falling and being impacted, as it may cause the user and/or the equipment.

Do not use a circular object or a material that is not the main one or accessories, clip using the device immediately, as it may cause damage to the device.

1. Do not use the device and wear accessories in an unstable situation, and prevent them from falling and being impacted, as it may cause the user and/or the equipment.

2. Do not use a circular object or a material that is not the main one or accessories, clip using the device immediately, as it may cause damage to the device.

3. Do not use the device for purposes other than those described in the user manual. This can cause the user and/or the equipment.

4. Do not disassemble, repair, or modify the device. This can cause harm, electric shock, and electric damage as well as increasing the possibility of purchasing a fake or an imitation device.

5. This is a normal blood glucose sensor; the device is powered on. There is no need to worry.

After using this product, temporary redness may occur due to the irritation. When changing the sensor, you have to wait with the adhesive tape (only L&L) is wet and probably to irritate. The irritation should last after the removal of the sensor. However, if the symptoms appear, please consult a professional doctor.

Common Questions

Q Can I use the outside of the cup?


A Yes, it can be used on the face and body. The available areas for use vary slightly depending on the attachment used. Please refer to "Notes and Functions" (p.17) for more detailed information.

Q What is the approximate usage time?

A Operate for 10-15 minutes on same area, and use the device for a total of approximately 20 minutes per day. Although it can be used every day, please rest/recharge the device on the same area.


Q Can I start with hot conditions or massage my cause skin problem. Therefore, please use the device on skin areas.

A Please use the device with hot conditions or massage may cause skin problems. Therefore, please use the device on skin areas.




No open flame
(cigarette, alcohol, etc.)

- Do not use the device if the AC adapter is broken when connected to an outlet, and it may cause a fire or electric shock.

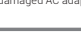

No charging
(charging station)


- Do not attach gel pads or gels to the charging port or charging terminal.
- Do not charge the device on the bed/room under any circumstances.
- Do not use the device when the plug is not on the same floor as the bed, mattress, or in a bath or shower.
- Do not use in a shower.



No high temperature

- Please operate in a room where the device is neither too hot, too cold, nor too dry. Use at room temperature, and do not place it near a fire, heat, or high temperature, and do not place it near a heater.
- Please be follow these precautions may result in electric shock, fire, and explosion.

Prohibit



Do not use power sockets.

When directly plugging the AC adapter into a wall socket, do not use metal contacts, as it may cause a fire, electric shock, or malfunction.

9

Weak or No AC stimulation

When the level button is pushed to increase the stimulation level, "Weak" or "No stimulation" appears on the LCD screen and the electrode pads and electrode pins are in close contact with the skin, not the skin. Be sure to apply enough gel when using the electrode pads.

Dist stimulation is too intense and painful

When the level button is pushed to increase the stimulation level, "Dist stimulation is too intense and painful" appears on the LCD screen. This may happen when the electrode pads and electrode pins are in close contact with the skin, not the skin. Be sure to apply enough gel when using the electrode pads.

When the intensity of DMS stimulation varies depending on the individual or the application area, adjust the level with a reasonable range. If the user experiences any abnormal sensations, stop using the device immediately and consult a doctor.

WARNING

! **Always follow these instructions:**

- Make sure that the main unit and charging terminal of the charger do not come into contact with water to prevent electric shock and fire.
- Use only the charging cable, AC adapter, accessories, and attachments that come with the power tool and the product.
- Use AC power of 100 VAC/1. If using the device continuously, be sure to confirm the shape of the plug and voltage before use to prevent fire, electric shock, or malfunctions.
- Before use, check whether there are any abnormalities in the pin and sockets. If you experience pain, pain abnormalities, or injuries during use, stop using the device immediately and consult a doctor to prevent personal accidents.

Quality and Dimensions

Main Unit	Approximately 129g	Approximately 76x40x108mm
Brush type attachment	Approximately 84g	Approximately 80x100x100mm
Thin type attachment	Approximately 84g	Approximately 100x100x100mm
Body attachment	Approximately 74g	Approximately 80x100x100mm

Materials

Main Unit: ABS resin (chrome plated), stainless steel, copper

Brush type attachment: ABS resin, chrome, stainless steel, copper, polycarbonate



Thin type attachment: ABS resin, chrome, stainless steel, copper, polycarbonate

Body attachment: ABS resin (chrome plated), PPS coating, stainless steel, copper, polycarbonate

Charging dock: ABS resin, slip, copper, polycarbonate

Power Consumption	Approximately 5W (when charging)	Lowest level: 1W (at least for a short period of time and you will get a warning. Adapt the level to your study programme)
Battery	Built-in lithium-ion battery 7.3V1700mAh	
Time	15-minute	
Battery Life⁺	Up to approximately 11 hours (when using level 1) Approximately 7 hours (when using level 5)	<p>⚠ If you have a laptop operation on the skin surface, reduce the level and move the device far from the skin surface. The local skin cooling system, using the gel to cool down the laptop and the skin.</p>

Note that the battery life may vary depending on the actual usage and environment of the device.
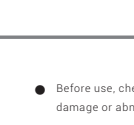
<div data-bbox="213 3038 263 3091">  </div> <div data-bbox="266 3038 441 3091"> <h2>WARNING</h2> </div> <div data-bbox="213 3091 441 3178">  <ul style="list-style-type: none"> Before each use, check the device for any damage, plug any part of the accessories to any damage or abnormalities to avoid any problems and damage. Remove any metal or inaccurate objects (electronic, springs, etc.) before using the device to avoid accidents and injuries. When using the device on the skin, do not touch or maintain prolonged contact (more than 10 min) and/or causing any damage. If you are using medical treatment, avoid using the device near the affected area and all of other medical devices while wearing, plug, or use the device. Do not bring children or animals and/or pets near it. (This area is not affected by the electrical field of the device.) </div>	<div data-bbox="445 3038 674 3091"> <h2>Storage methods</h2> </div> <div data-bbox="445 3091 674 3178"> <ul style="list-style-type: none"> ● If after use, please wipe the device and attachment dry and store them in a dry and well-ventilated area on the charging dock. Please note that the charging dock is not waterproof, and placing the device in a wet state on the charging dock may cause malfunctions. Using a cotton ball to clean the electrical contacts of the device and charging dock to remove any water or dust. ● The device is not used for a long time, it should be cleaned and maintained before storage. ● Turn off the device and unplug the AC adapter from the wall socket when not in use. ● Please do not store the device in high-temperature and high-humidity areas such as near open flames, in the bathroom, or in a hot car. ● Store the device in a clean, cool, and well-ventilated area out of reach of children. ● The device is not used for a long time, please use it fully each time. </div>
--	---

[illegible]

Prohibit

- Do not use a generic touch screen that does not apply this device properly. The touch screen and electrode must be used as directed.
- Do not use the device on any of the following individuals (just use the device):
 - Those with skin sensitivity. Those who do have sensitive skin, band-aids, and those who wear nail or types of jewelry (just use an electrode tape or an electrode).
- Do not use the device continuously for a long time. The device may become hot with extended use. Heat for an extended use has been shown to damage organs.

→ Please note the device is the best for the use of the device.

⚠️ The device may be used on the device.

15

40

Before using the device, please confirm the following:


If any of the accessories are missing from the package, please contact customer service.

The package includes:

- 1 device,
- 3 attachments (brush-type, for face, for body),
- 1 charging dock,
- 1 AC adapter,
- 1 instruction manual.

④ The design of the device and accessories may change without prior notice to improve quality.

● For spot care, apply the device on the black spots for 3-5 seconds each, following the sequence below. Repeat the entire sequence for 2-4 minutes.



Press down the device directly on the spot for 3-5 seconds. Repeat the entire sequence for 2-4 minutes.

④ When using the device on the face, please do not press too hard on the blemished area, cheeks, or chin, as they may cause discomfort. Please use the device on each part separately. When using the device on the head/neckline, please make sure that it may cause a strong stimulation.

[illegible]

■ STEP1

For the front and top of the head, place the device in a zigzag pattern. (Usually move from one point to another.)

■ STEP2

For the sides of the head, start at the hairline and move towards the back of the head. (Usually move from one point to another.)

■ STEP3

For the back of the head, move the device from the upper side towards the collarbone and back down. (Usually move from one point to another.)

■ Charging dock

Charging terminal

Dedicated AC charger power (not sold sep.)

■ Dedicated AC charger

AC adapter

Micro-connector

⚠Charging dock and dedicated AC charger are not waterproof. Do not place them in wet conditions, and do not store them while wet. Absolutely do not charge them in the bathroom.

Use the Lifting Mode

Get Required

Use the Relaxing Mode

Get Required

Use with the brush-shaped electrode for scalp and skin care.

Use with the brush-shaped electrode with a lifting sensation.

The skin lifting pain stimulates the skin and muscle under the skin, and improves the skin elasticity. Care should be on the scalp, face, and forehead.

Use with the facial muscle electrode and body care.

The long DMS pads contract the muscles, care for the facial muscle and body care. Care should be on the face and body (especially back) with various accessories.


Before feel a strong stimulation, start at a low level (level 1) and use the device for a short period of time.

Although the device can be used on damp scalp and skin with toner applied, please use it at a moderate stimulation level.

Please do not use a device on a child who does not fully understand the instructions. If you have any questions, please ask your doctor or contact us before using the device. When you are driving the car, please do not use the device to ensure safe driving. Please do not use the device on the nose, neck, temples, or the area around the eyes (the area marked as X).

21

34



When using the face or body attachments, the sensation of electrical stimulation may vary depending on the area of the body and the intensity of the treatment. When using the face attachment, the intensity of the treatment should be adjusted to the individual's skin type. When using the body attachment, the intensity of the treatment should be adjusted to the individual's body type.

When using the face attachment, the device should be held against the face with the arrow pointing towards the face. When using the body attachment, the device should be held against the body with the arrow pointing towards the body.


To achieve more effective facial care, you can combine two different attachments for different purposes.

For achieving more effective facial care, you can combine two different attachments for different purposes. For the face attachment, the device should be held against the face with the arrow pointing towards the face. For the body attachment, the device should be held against the body with the arrow pointing towards the body.

Attachment	Intensity	Duration
Face, head, shoulders & neck	10-15 (mmHg)	10-15 (minutes)
Face, head, shoulders & neck	5 (mmHg)	5 (minutes)


When using the face attachment, the device should be held against the face with the arrow pointing towards the face. When using the body attachment, the device should be held against the body with the arrow pointing towards the body.

- 6 After completing the one press and the power button for approximately 3 seconds is turn off the device.
 - The water will not automatically leak off. It will remain wet.
 - Please use our daily storage product after the battery.
- 7 Please note the following regarding the waterproof function.

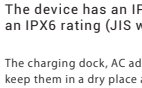


The device has an IPX5 waterproof rating, and the attachments have an IPX5-rated (30 minute protection level).

The charging dock, IPX5 cables, and attachments do not have a waterproof rating. Please keep them in a dry place and do not use them in wet environments, including during charging in the rain.

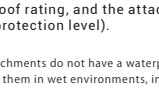

When the battery level is low, LED 1 will flash quickly 2 times every 3 seconds.




Charging time is approximately 2 hours, and operating time is approximately 7 hours when used as usual (on the basis that the attached, as a noise reduction of 10 dB).
 (Detail) About charging and playing time, see the following on the battery power and charging environment.


- 2 Select an accessory and attach it to the main unit. The accessory is magnetic. Insert it by aligning the mounting part of the accessory with the magnet part of the main unit (see diagram below).
 - The device will automatically switch to the appropriate mode (full HD or deep HD) mode based on the accessory being used.
 - To remove the accessory, hold it from the side and pull it off.
- 3 Connect the AC adapter in the charging dock then place the device into the charging dock with the charging terminals connected.
 - Before charging, wear your water droplets off the device and charging terminals. Moisture on the charging terminals can cause malfunctions.
- 4 When charging is complete, unplug the AC adapter from the wall outlet.

- 2 Turn on the device and wait at least 10 minutes for your eyes to adapt to the light. Turn the device on until it shows "TEST 1" on the display.
- 3 If you experience any abnormal reactions such as red, itchy, or pain during the test, immediately stop using the device and seek medical attention.
- 4 Wait for 24 hours after the test and check for any abnormalities such as pain or heat sensation.
- 5 Obtain the temporary readings may occur on the skin after using the device due to the stimulation. The above test results may be different if you wear glasses or contact lenses, please consult a doctor.

Please measure any abnormal reaction during the patch test using the device immediately.



- 1 Use the device in the same way for 5-15 minutes, and do not use it for more than 30 minutes per day. If you experience any abnormal reactions resulting from extended use, stop using the device immediately.
- 2 The device will warn you by emitting a warning tone and the warning tone will be shown. When the device will warn you by emitting a warning tone, stop using the device immediately.
- 3 When the device does not emit a warning tone, the same level of stimulation can be used for the next use after 24 hours.
- 4 Avoid excessive use and gradually lower the device voltage if you use the device frequently in any particular skin condition or skin type.
- 5 If you experience any discomfort or abnormal sensation on your scalp or skin, stop using the device immediately and consult a doctor.

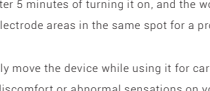

- 1 Referencing the device, please ensure that your scalp and skin are clean and free of cosmetics, sunscreen, hair cream, etc. Before use, please check whether there are growths or inflammation on the scalp and skin. Do not use the device on areas with abnormalities.